

Drainage Bag Closed System

Urimaax Drainage Bag Closed System is a 2000ml drainage bag with a bottom T-bar drain with large finger control and a universal stepped catheter connector. This system has a needle free sampling point and hanger holes that fit universal plastic or metal hangers. Ideal for use with Urimaax Drainage Bag Accessories.



One-Handed Application

Large finger control allows for one-handed use




Easy to Read

100ml graduated measurements for easy identification



Fast and Safe Sampling

Needle free sampling point allows for urine sampling while maintaining a closed system

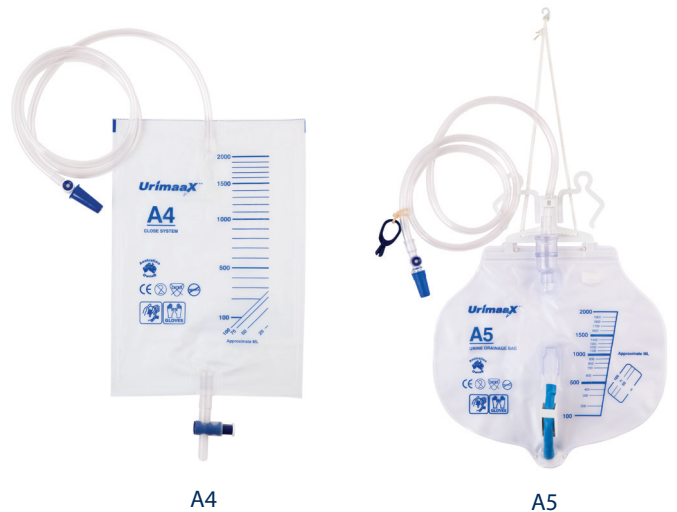

A4

A5

- Suitable for:**
- Suprapubic or indwelling catheters
 - One-handed application
 - Most bed shapes

REF	Type	Capacity	Sterility	Pcs/Unit	Units/Ctn
10013104	A4	2000ml	Sterile	50	1
10013105	A5	2000ml	Sterile	20	1

At Sentry, we are dedicated to upholding sustainable principles to safeguard our planet, people, and future, and we strive to consistently act responsibly in doing so.

